

Multiple shooting methods



Handheld shooting



Remote control with a tripod and extended exposure hand brake line

Technical Parameter **Portable X-ray machine**

| | | |
|-----------------------------|----------------------------|---------------------------|
| Tube Voltage 60KV | Highest Voltage 70KV | Tube Current 2mA |
| Exposure Time 0.01-2s | Focus 0.4mm | Anode Target Angle 12° |
| Total Filtration 1.0mmAl | Battery Capacity 950mAh | |

**Portable X-ray**

A smart dental digital solution

داینا دنت
نماینده انحصاری فیوژن در ایران
fussengroup.ir


Fussen Technology Co.,Ltd

1F,3#Building ,Tongfuyu Industrial Zone, Taoyuan Street,Nanshan District, Shenzhen,China

Tel:+86 755 2167 0890 - 813


Website: www.fussengroup.com

Email:sales@fussenct.com

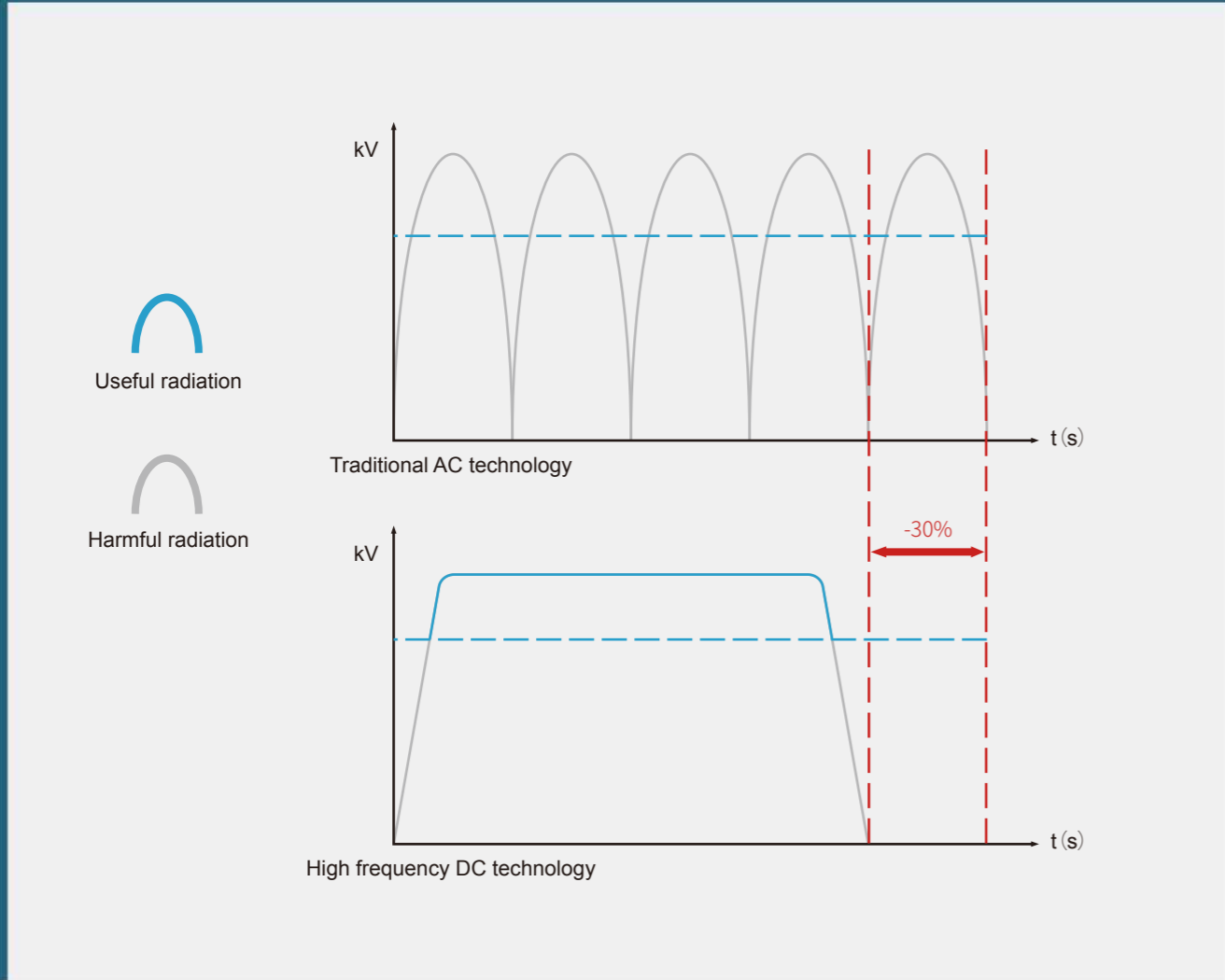

High frequency DC technology

Ultra-low radiation, the exposure dose is reduced by more than 30% compared with the standard generator, maximizing effective rays and greatly reducing the radiation intensity

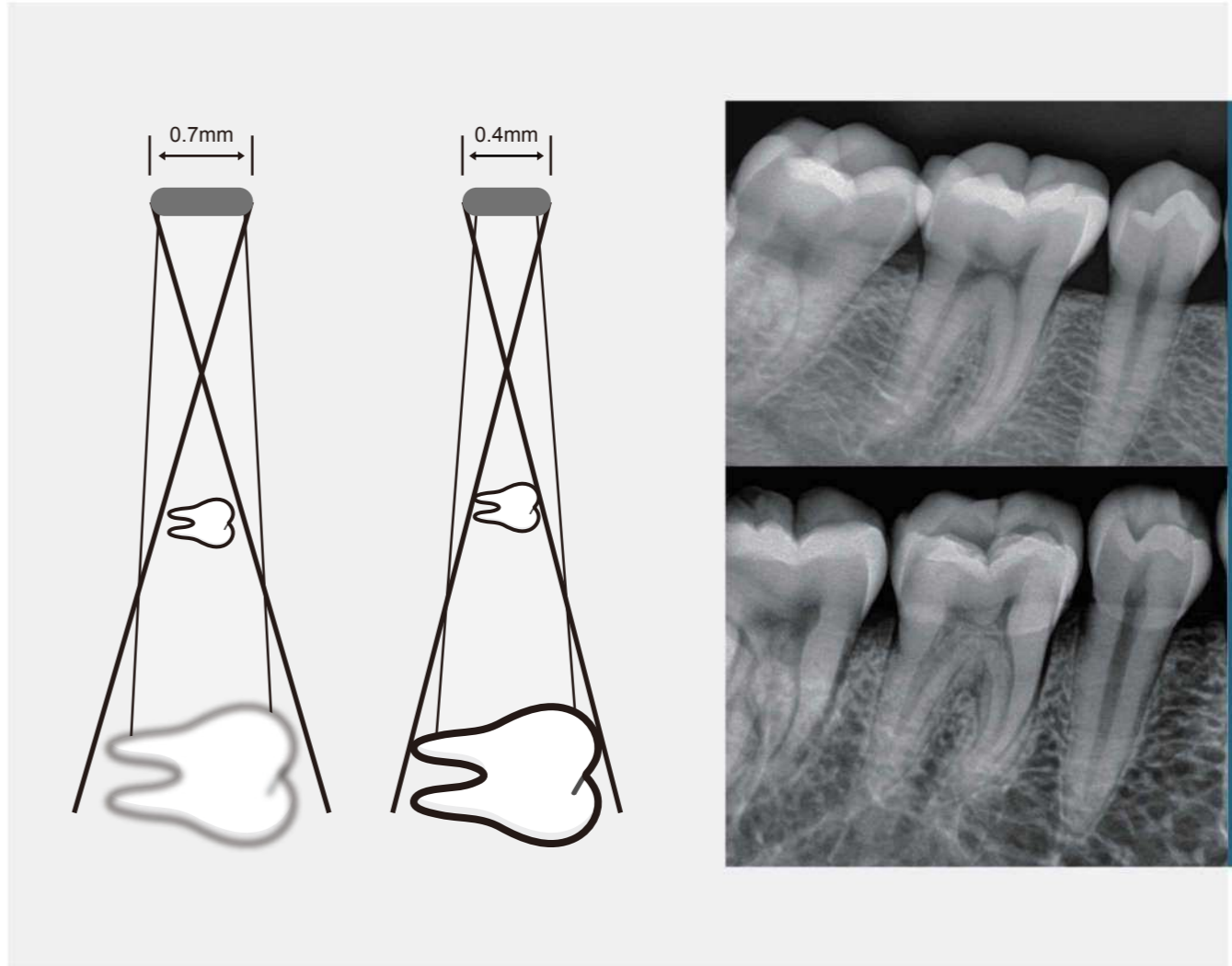

Small focus to get high definition images

High-voltage generator with 0.4mm focus, high-definition imaging, enjoy top quality



High capacity lithium battery

950 mAh ultra high capacity battery, battery can be replaced, low cost of use



Humanized and easy operation

Shoot different teeth, fine-tune the shooting time, and achieve personalized and precise control

